

## Release Notes

**mint.jar 1.24.08**  
(Released 13/12/2022)

### DICOM SR to Encapsulated PDF Node


Added in a new node to convert a DICOM SR to an encapsulated PDF DICOM file. The headings in the report are customisable. The original SR passes through the node and the created file is sent to the defer channel.

### PDF Output Node

Added in a new node to extract the PDF file from an encapsulated PDF DICOM file and write the PDF to the specified location with a .pdf extension.

### DICOM File Writer Node

Updated the DICOM File Writer node to allow a custom flat folder name definition. The sections of the folder (separated by \_) are defined as a tag list.

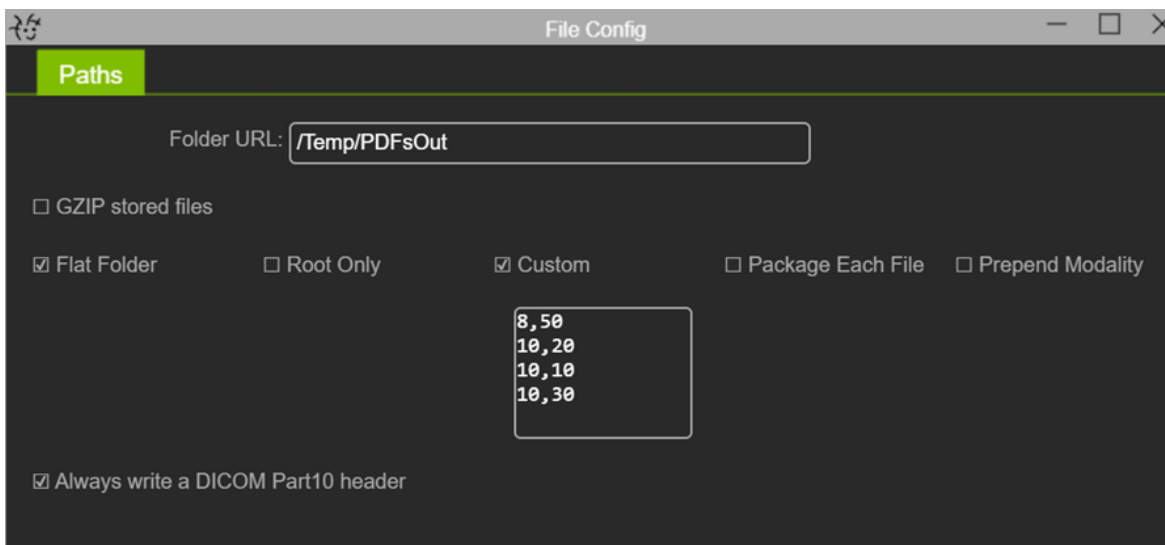

For example, the config above would include the accession number with patient id, name and dob in the folder name as {root}/study/\$ACC\_NUM\_\$PID\_\$PAT\_NAME\_\$PAT\_DOB